

"The Clinical Advantage"™

NEW

Announcing the birth of the next generation
Ultrasound & Echocardiography Tables



BIODEX

www.biodes.com

1-800-224-6339

In New York and Int'l, call 631-924-9000

Accommodates
BARIATRIC
Patients



BIODEX MEDICAL

NEW

Sound Pro Combination Table

Dual-purpose table accommodates ultrasound and echocardiography procedures

- Ergonomically designed to help prevent musculoskeletal injuries
- Extra wide top with 500-lb patient capacity – perfect for bariatric patients
- Motor-controlled adjustments – including Fowler back, height, Trendelenburg motions and auto level
- Cardiac scanning cutout features a patented swing-out cushion that doubles as an armrest
- Adjustable sonographer cutout allows comfortable scanning positioning
- Vascular scanning arm board can be mounted to either side of the table
- Side rails pivot out of the way to allow unobstructed patient access
- Fully equipped with table accessories including hand and foot controllers, pivoting side rails, vascular scanning arm board, retractable stirrups, IV pole and paper dispenser



The Sound Pro Combination Table is a versatile table solution, combining the imaging features of both the Ultra Pro and Echo Pro, making it ideal for OB/GYN, general ultrasound and echocardiography. Fully equipped with ultrasound and echocardiography capabilities, the Sound Pro is extra wide with

500-lb patient capacity and a vast height range down to 23", easily accommodating wheelchair transfers, patients with ambulatory difficulties and the bariatric population.

Ergonomically designed to suit both patient and sonographer, the table's height range allows seated or standing work positions while the pivoting side rails and adjustable sonographer cutout cushion allow comfortable shoulder, arm and hand positions when scanning. There are no barriers between sonographer and patient, allowing correct scanning posture to prevent musculoskeletal injury. The cardiac scanning cutout allows open access to the patient's left thorax area for an unobstructed apical approach. The motorized Fowler back is infinitely adjustable up to 80° via hand or foot controller. The table is so accommodating, stress-echo patients are off the treadmill, onto the table and into position for imaging quickly and comfortably.

For the patient, comfort is assured through lengthy procedures with cushioned arm supports, an amply thick mattress and a multitude of comfort-driven design details. For example, the unique swing-out cushion of the cardiac scanning cutout doubles as an armrest. The leg section drops down to 40° for patients with circulatory issues. The extra-wide tabletop comfortably accommodates the bariatric patient. In fact, the table is so accommodating, it reduces the total amount of time required to achieve a quality image.

The Sound Pro Combination Table is fully equipped... no additional options or accessories are required. Features include: 2"- and 3"-thick table padding covered in easy-to-clean Naugahyde®, one patient restraining strap with Velcro® hook and loop fasteners, five-inch swivel casters and a central locking system, accessible from either side of the table, to secure table position.

SPECIFICATIONS:

Dimensions: 71" l x 30" w (180.3 x 76.2 cm); 35.50" w (90.2 cm) with side rails

Drop-Down Leg Section: 12.6" l (32 cm)

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Tabletop: Three primary sections: torso section features sonographer and patient cutouts with Fowler back positioning; center section remains fixed and leg section drops down to 40° and 80° for stirrup access

Cardiac Scanning Cutout: with pivoting cushion that doubles as an armrest

Adjustable Sonographer Cutout: drops down or folds up for maximum access to the patient

Mattress: Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions, Fowler back and auto level motions

Foot Control: Activates height, Trendelenburg motions and Fowler back

Wheels: 5" (12.7 cm) swivel casters; central floor locking system

Patient Restraints: One body strap

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.


Weight: 350 lb (159 kg)

Power: 115 VAC or 230 VAC

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2

No. 601.1-M90 and EN 60601-1 standards.

Warranty: One year

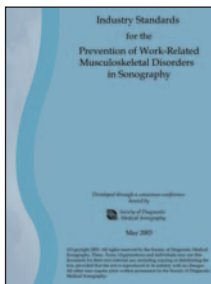

More than 89% of sonographers are scanning in pain and 20% of these professionals eventually experience a career-ending injury. On average, within five years of entering the profession, sonographers experience pain while scanning.

For more information on SDMS Industry Standards for the Prevention of Work-Related Musculoskeletal Disorders in Sonography, visit our

website at www.biomed.com

058-710 Table, Combination, Sound Pro, 115 VAC

Includes: hand and foot controllers, pivoting side rails, vascular scanning arm board, retractable stirrups, IV pole and paper dispenser.

058-715 Table, Combination, Sound Pro, 230 VAC

Includes: hand and foot controllers, pivoting side rails, vascular scanning arm board, retractable stirrups, IV pole and paper dispenser.

Replacement:

056-612 Table Paper, Hygenic, 12 rolls

18" w x 225' l (46 cm x 69 m)

BIODEX

Accommodates
BARIATRIC
Patients



BIODEX MEDICAL

NEW Echo Pro Echocardiography Table

Efficient positioning reduces sonographer injuries and shortens exam time.

- Ergonomically designed to help prevent musculoskeletal injuries
- Extra wide top with 500-lb patient capacity – perfect for bariatric patients
- Motor-controlled adjustments – including Fowler back, height, Trendelenburg motions and auto level
- Cardiac scanning cutout features a patented swing-out cushion that doubles as an armrest
- Adjustable sonographer cutout allows comfortable scanning positioning
- Vascular scanning arm board (optional) can be mounted to either side of the table
- Side rails (optional) pivot out of the way to allow unobstructed patient access
- Equipped with hand controller to actuate table positions
- Optional table accessories include foot controller, vascular scanning arm board, pivoting side rails, IV pole and paper dispenser



The Echo Pro Echocardiography Table is a dedicated solution for all echocardiography procedures. Extra wide with 500-lb patient capacity and a vast height range down to 23", this table will easily accommodate wheelchair transfers, patients with ambulatory difficulties and the bariatric population.

Ergonomically designed to suit both patient and sonographer, the table's height range allows seated or standing work positions while the pivoting side rails and adjustable sonographer cutout cushion allow comfortable shoulder, arm and hand positions when scanning. There are no barriers between sonographer and patient, allowing correct scanning posture to prevent musculoskeletal injury. The cardiac scanning cutout allows open access to the patient's left thorax area for an unobstructed apical approach. The motorized Fowler back is infinitely adjustable up to 80° via hand or foot controller. The table is so accommodating, stress-echo patients are off the treadmill, onto the table and into position for imaging quickly and comfortably.

For the patient, comfort is assured through lengthy procedures with cushioned arm supports, an amply thick mattress and a multitude of comfort-driven design details. For example, the unique swing-out cushion of the cardiac scanning cutout doubles as an armrest. The extra-wide tabletop comfortably accommodates the bariatric patient. In fact, the table is so accommodating, it reduces the total amount of time required to achieve a quality image.

The Echo Pro Echocardiography Table is equipped with a motor-controlled hand controller to actuate table positions, 2"- and 3"-thick table padding covered in easy-to-clean Naugahyde®, two patient restraining straps with Velcro® hook and loop fasteners, five-inch swivel casters and a central locking system, accessible from either side of the table, to secure table position.

Additional accessories are available, including foot controller, vascular scanning arm board, pivoting side rails, IV pole and paper dispenser.

SPECIFICATIONS:

Dimensions: 73" l x 30" w (185.4 x 76.2 cm); 35.5" w (90.2 cm) with optional side rails

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Tabletop: Two primary sections: torso section: features sonographer and patient cutouts with Fowler back positioning, and lower body section remains fixed

Cardiac Scanning Cutout: with pivoting cushion that doubles as an armrest

Adjustable Sonographer Cutout: drops down or folds up for maximum access to the patient

Mattress: Torso section 2" (5 cm) thick, lower body section 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions, Fowler back and auto level motions

Foot Control (optional): Activates height, Trendelenburg motions, and Fowler back

Wheels: 5" (12.7 cm) swivel casters; central floor locking system

Patient Restraints: Two body straps

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.


Weight: 335 lb (152 kg)

Power: 115 VAC or 230 VAC

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2

No. 601.1-M90 and EN60601-1 standards.

Warranty: One year

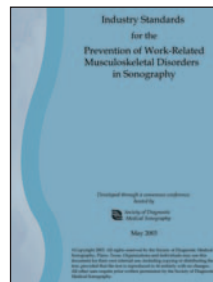

More than 89% of sonographers are scanning in pain and 20% of these professionals eventually experience a career-ending injury. On average, within five years of entering the profession, sonographers experience pain while scanning.

For more information on SDMS Industry Standards for the Prevention of Work-Related Musculoskeletal Disorders in Sonography, visit our

website at www.biomed.com

058-700 Table, Echocardiography, Echo Pro, 115 VAC
Includes Hand Controller.

*All options available separately**

058-705 Table, Echocardiography, Echo Pro, 230 VAC
Includes Hand Controller.

*All options available separately**

Optional:

058-740 Foot Controller

058-631 Side Rails, Pivoting

058-735 Arm Board, Vascular Scanning

058-656 IV Pole

058-610 Dispenser, Paper
with Cutter Strap

Replacement:

056-612 Table Paper, Hygienic, 12 rolls

BIODEX



BIODEX MEDICAL

NEW

Ultra Pro Ultrasound Table

Sonographer-friendly design.

- Ergonomically designed to help prevent musculoskeletal injuries
- Contoured tabletop design for optimal probe positioning
- 500-lb patient capacity
- Open frame design for greater accessibility
- Motor-controlled adjustments – including Fowler back, height, Trendelenburg motions and auto level
- Vascular scanning arm board (optional) can be mounted to either side of the table
- Side rails (optional) pivot out of the way to allow unobstructed patient access
- Equipped with hand controller to actuate table positions
- Optional table accessories include foot controller, retractable stirrups, vascular scanning arm board, pivoting side rails, IV pole and paper dispenser



The Ultra Pro Ultrasound Table is a dedicated table for all general ultrasound procedures including OB/GYN. Choose from an assortment of table options to customize the table to suit your department needs.

Contoured tabletop design, with 500-lb patient capacity, maximizes patient scanning. With a vast height range down to 23", the Ultra

Pro easily accommodates wheelchair transfers and patients with ambulatory difficulties.

Ergonomically designed to suit both patient and sonographer, the table's height range allows seated or standing work positions while the pivoting side rails allow comfortable shoulder, arm and hand positions when scanning. There are no barriers between sonographer and patient, allowing correct scanning posture to prevent musculoskeletal injury. The open frame design allows sonographers to get closer to the patient, especially useful for carotid imaging or long leg arterial scanning. The motorized Fowler back is infinitely adjustable up to 80° via hand or foot controller.

For the patient, comfort is assured through lengthy procedures with an amply thick mattress and a multitude of comfort-driven design details. For example, the retractable stirrups can extend extra long to accommodate taller or larger patients; the leg section drops down to 40° for patients with circulatory issues and 80° for stirrup access. Biodex tables are so accommodating, they actually reduce the total amount of time required to achieve a quality image.

The Ultra Pro Ultrasound Table is equipped with a motor-controlled hand controller to actuate table positions, 2"- and 3"-thick table padding covered in easy-to-clean Naugahyde®, two patient restraining straps with Velcro® hook and loop fasteners, five-inch swivel casters and a central locking system, accessible from either side of the table, to secure table position.

Additional accessories are available, including foot controller, retractable stirrups, vascular scanning arm board, pivoting side rails and IV pole.

SPECIFICATIONS:

Dimensions: 70" l x 27" w (185.4 x 68.6 cm); 32.5" w (82.5 cm) with side rails; 18" w (45.72 cm) at narrow section

Drop-Down Leg Section: 12.6" l (32 cm)

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Tabletop: Three primary sections: torso section features contoured design with Fowler back positioning, center section remains fixed and leg section drops down to 40° and 80° for stirrup access

Mattress: Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions, Fowler back and auto level motions

Foot Control (optional): Activates height, Trendelenburg motions, and Fowler back

Wheels: 5" (12.7 cm) swivel casters; central floor locking system

Patient Restraints: Two body straps

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.


Weight: 335 lb (152 kg)

Power: 115 VAC or 230 VAC

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2

No. 601.1-M90 and EN 60601-1 standards.

Warranty: One year

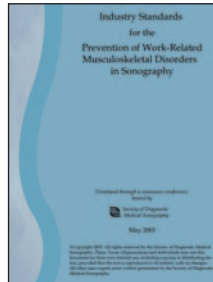

More than 89% of sonographers are scanning in pain and 20% of these professionals eventually experience a career-ending injury. On average, within five years of entering the profession, sonographers experience pain while scanning.

For more information on SDMS Industry Standards for the Prevention of Work-Related Musculoskeletal Disorders in Sonography, visit our

website at www.biodex.com

058-720 Table, Ultrasound, Ultra Pro, 115 VAC

Includes Hand Controller.

*All options available separately**

058-725 Table, Ultrasound, Ultra Pro, 230 VAC

Includes Hand Controller.

*All options available separately**

Optional:

058-740 Foot Controller

058-631 Side Rails, Pivoting

058-651 Stirrups, Retractable

058-735 Arm Board, Vascular Scanning

058-656 IV Pole

058-610 Dispenser, Paper

with Cutter Strap

Replacement:

056-612 Table Paper, Hygenic, 12 rolls

18" w x 225' l (46 cm x 69 m)

BIODEX



BIODEX MEDICAL

Deluxe Ultrasound Table

Designed for all Ultrasound and OB/GYN procedures.

- Ergonomically designed to help prevent musculoskeletal injuries
- Contoured tabletop design for optimal probe positioning
- Open frame design for greater accessibility
- Manual positioning for Fowler back, height and Trendelenburg motions
- Side rails (optional) fold down, out of the way to allow unobstructed patient access



The Deluxe Ultrasound Table features three primary positions: horizontal, chair and Trendelenburg.

A convenient and easy-to-operate hydraulic foot pump adjusts table height during patient transfers and scanning procedures. Trendelenburg positioning is easily performed with a smooth-turning hand crank. Ankle straps, located at the underside of the leg section, prevent patient sliding during Trendelenburg positioning.

Fowler positioning is accommodated with preset stops at 0°, 15°, 35°, 55° and 80°. The tabletop's hourglass shape narrows to 19 inches, providing comfortable scanning probe positioning at and below the plane of the tabletop. The leg section unlocks and swings down for chair positioning, or for positioning with optional stirrups.

Patient comfort during lengthy procedures is assured by two-inch thick table padding covered in durable and easy-to-clean Naugahyde®. Two adjustable straps with Velcro® hook and loop fastening are provided for increased patient stability and comfort. Five-inch locking swivel casters provide smooth and secure table positioning. Optional folding side rails drop down and out of the way for unobstructed access to the patient.

SPECIFICATIONS:

Dimensions: 75" l x 24" w (191 x 61 cm); 30" w (76 cm) with side rails;

19" w (48 cm) at narrow section

Drop-Down Leg Section: 14.5" l (37 cm), drops down for stirrup access

Height Adjustable: 29.5" to 37" (75 to 94 cm)

Tabletop: Three primary sections; torso section features contoured design with Fowler back positioning, center section remains fixed and leg section drops down for stirrup access

Mattress: 2" (5 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 29.5" to 37" (75 to 94 cm), hydraulic foot pump

Trendelenburg: 0° to 15°, hand crank

Fowler Back: Manually adjusted in five preset increments, 0°, 15°, 35°, 55° and 80°

Wheels: 5" (13 cm) locking casters

Patient Restraints: Two body straps

Finish: Powder coat

Patient Capacity: 350 lb (159 kg); weight tested to four times the patient load rating.

Weight: 335 lb (152 kg)

Warranty: One year

056-613 Table, Ultrasound, Deluxe
with Convenience Package

Includes: Side Rails, Stirrups and IV Pole

056-609 Table, Ultrasound, Deluxe

*Includes: Side Rails**

056-605 Table, Ultrasound, Deluxe

*All options available separately.**

Optional:

056-608 Stirrups

056-607 IV Pole

058-610 Dispenser, Paper
with Cutter Strap

Replacement:

056-612 Table Paper, Hygienic, 12 rolls
18" w x 225' l (46 cm x 69 m)

**Note: Options are installed by the customer.*

Compare the **NEW** Ultrasound tables from Biodex

Which table is right for you?

NEW

Sound Pro Combination Table



SPECIFICATIONS:

Dimensions: 71" l x 30" w (180.3 x 76.2 cm); 35.50" w (90.2 cm) with side rails
Drop-Down Leg Section: 12.6" l (32 cm)
Height Adjustable: 23" to 39" (58.4 to 99 cm)
Tabletop: Three primary sections: torso section features sonographer and patient cutouts with Fowler back positioning; center section remains fixed and leg section drops down to 40° and 80° for stirrup access

Cardiac Scanning Cutout: with pivoting cushion that doubles as an armrest

Adjustable Sonographer Cutout: drops down or folds up for maximum access to the patient

Mattress: Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions,

Fowler back and auto level motions

Foot Control: Activates height, Trendelenburg motions

and Fowler back

Wheels: 5" (12.7 cm) swivel casters; central floor locking system

Patient Restraints: One body strap

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.

Weight: 350 lb (159 kg)

Power: 115 VAC or 230 VAC

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2

No. 601.1-M90 and EN 60601-1 standards.

Warranty: One year

058-710 Table, Combination, Sound Pro, 115 VAC
Includes: hand and foot controllers, pivoting side rails, vascular scanning arm board, retractable stirrups, IV pole and paper dispenser.

058-715 Table, Combination, Sound Pro, 230 VAC
Includes: hand and foot controllers, pivoting side rails, vascular scanning arm board, retractable stirrups, IV pole and paper dispenser.

Replacement:

056-612 Table Paper, Hygenic, 12 rolls
 18" w x 225' l (46 cm x 69 m)

NEW

Echo Pro Echocardiography Table



SPECIFICATIONS:

Dimensions: 73" l x 30" w (185.4 x 76.2 cm); 35.5" w (90.2 cm) with optional side rails
Height Adjustable: 23" to 39" (58.4 to 99 cm)
Tabletop: Two primary sections: torso section: features sonographer and patient cutouts with Fowler back positioning, and lower body section remains fixed

Cardiac Scanning Cutout: with pivoting cushion that doubles as an armrest

Adjustable Sonographer Cutout: drops down or folds up for maximum access to the patient

Mattress: Torso section 2" (5 cm) thick, lower body section 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions,

Fowler back and auto level motions

Foot Control (optional): Activates height, Trendelenburg motions, and Fowler back

Wheels: 5" (12.7 cm) swivel casters; central floor locking system

Patient Restraints: Two body straps

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.

Weight: 335 lb (152 kg)

Power: 115 VAC or 230 VAC

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2

No. 601.1-M90 and EN60601-1 standards.

Warranty: One year

058-700 Table, Echocardiography, Echo Pro, 115 VAC
Includes Hand Controller.

*All options available separately**

058-705 Table, Echocardiography, Echo Pro, 230 VAC
Includes Hand Controller.

*All options available separately**

Optional:

058-740 Foot Controller

058-631 Side Rails, Pivoting

058-735 Arm Board, Vascular Scanning

058-656 IV Pole

058-610 Dispenser, Paper

with Cutter Strap

Replacement:

056-612 Table Paper, Hygenic, 12 rolls
 18" w x 225' l (46 cm x 69 m)

**Note: Options are installed by the customer.*

NEW

Ultra Pro Ultrasound Table



SPECIFICATIONS:

Dimensions: 70" l x 27" w (185.4 x 68.6 cm); 32.5" w (82.5 cm) with side rails;

18" w (45.72 cm) at narrow section

Drop-Down Leg Section: 12.6" l (32 cm)

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Tabletop: Three primary sections: torso section features contoured design with Fowler back positioning, center section remains fixed and leg section drops down to 40° and 80° for stirrup access

Mattress: Torso section 2" (5 cm) thick, center and leg sections 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®

Upholstery Color: Graphite

Motions:

Height Adjustable: 23" to 39" (58.4 to 99 cm)

Trendelenburg: 0° to ±15°

Fowler Back: 0° to 80° infinitely adjustable

Controls:

Hand Control: Activates height, Trendelenburg motions,

Fowler back and auto level motions

Foot Control (optional): Activates height, Trendelenburg motions, and Fowler back

Wheels: 5" (12.7 cm) swivel casters; central floor locking system

Patient Restraints: Two body straps

Finish: Powder coat

Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.

Weight: 335 lb (152 kg)

Power: 115 VAC or 230 VAC

Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2

No. 601.1-M90 and EN 60601-1 standards.

Warranty: One year

058-720 Table, Ultrasound, Ultra Pro, 115 VAC
Includes Hand Controller.

*All options available separately**

058-725 Table, Ultrasound, Ultra Pro, 230 VAC
Includes Hand Controller.

*All options available separately**

Optional:

058-740 Foot Controller

058-631 Side Rails, Pivoting

058-651 Stirrups, Retractable

058-735 Arm Board, Vascular Scanning

058-656 IV Pole

058-610 Dispenser, Paper

with Cutter Strap


Replacement:

056-612 Table Paper, Hygenic, 12 rolls
 18" w x 225' l (46 cm x 69 m)

**Note: Options are installed by the customer.*

I can give you
Five Good Reasons
to buy a New Biodex Ultrasound Table...

1. Five-hundred pound patient capacity
2. Power Fowler back
3. Removable arm board
4. Auto-level
5. Superior ergonomic design

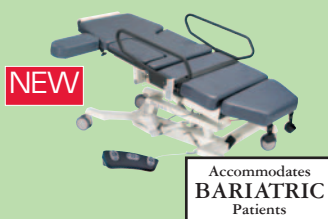

Sound Pro
Combination
Table

Combining all the benefits and features of the Ultra Pro and Echo Pro into one universally adaptable table. Perfect for the bariatric patient.



Echo Pro
Echocardiography
Table

Comfortable access and efficient positioning reduce sonographer injuries, shorten exam time and improve results



Ultra Pro
Ultrasound
Table

"Sonographer friendly" design helps prevent musculoskeletal disorders

BIODEX

Biodex Medical Systems, Inc.

20 Ramsay Road, Shirley, New York, 11967-4704, Tel: 800-224-6339 (In NY and Int'l, call 631-924-9000), Fax: 631-924-9241, Email: info@biodex.com, www.biodex.com



BIODEX

Biodex Medical Systems, Inc.