



A view of the normal inferior turbinate with a traditional halogen light source.

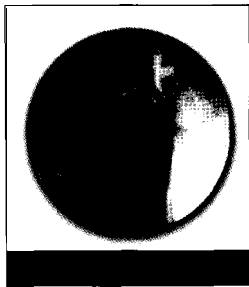


The normal inferior turbinate as illuminated by the new Welch Allyn Solarc light source.*

*Photos have been reproduced to simulate the results you can expect.

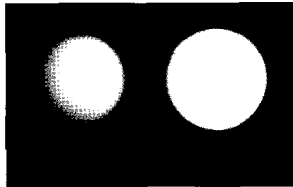

Welch Allyn's patented Solarc arc lamp provides white, xenon-quality light (5,500°K) for more accurate color rendition compared to halogen systems (3,200°K). The proprietary lamp design focuses more light into the fiberoptic illumination bundle for brighter light output. The Solarc has an average lamp life of 350 hours and, at 2.2 lbs., is completely portable.

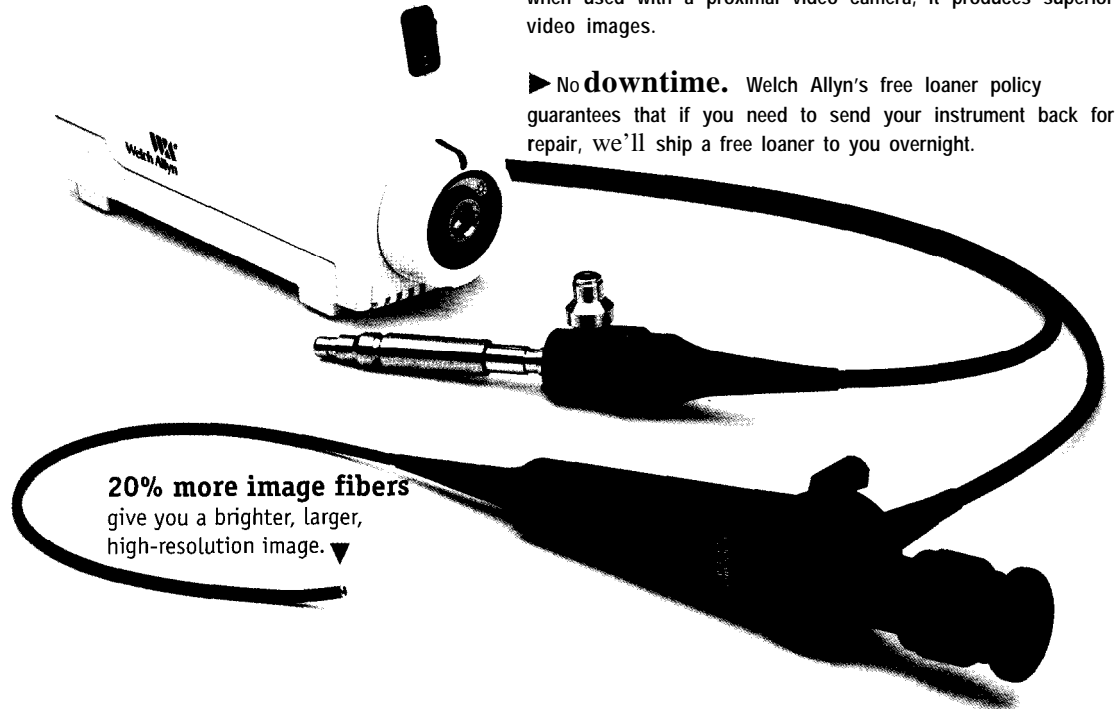

▶ At only 2.2 lbs., the Solarc has all the portability you need for postoperative exams in the hospital.

▶ The Solarc light source is also **compatible with other instrumentation**, including rigid sinusscopes, other rhinolaryngoscopes and fiberoptic headlights. Furthermore, when used with a proximal video camera, it produces superior video images.

▶ **No downtime.** Welch Allyn's free loaner policy guarantees that if you need to send your instrument back for repair, we'll ship a free loaner to you overnight.

20% more image fibers give you a brighter, larger, high-resolution image. ▼

SPECIFICATIONS

| | |
|--------------------------|-----------|
| ANGLE OF VIEW | 75° |
| DEPTH OF FIELD | 3-50 MM |
| DEFLECTION [UP/DOWN] | 130°/130° |
| BENDING SECTION DIAMETER | 3.4 MM |
| INSERTION TUBE DIAMETER | 3.5 MM |
| INSERTION TUBE LENGTH | 30 CM |
| IMMERSIBLE | YES |

| ORDERING INFORMATION | MODEL# |
|---|--------|
| RL-150 RHINOLARYNGOSCOPE AND Sdarc LIGHT SOURCE | 65053 |
| RL-150 RHINOLARYNGOSCOPE ONLY | 65050 |
| Sdarc LIGHT SOURCE ONLY | 46100 |
| STORAGE HANGER | 65003 |
| DISINFECTION TRAY | 65002 |

Specially designed disinfection tray makes reprocessing more convenient. ▶



◀ An optional storage hanger easily secures the instrument when not in use, **extending the life of the RL-150** by protecting it against damage and loss.

WELCH ALLYN 4341 STATE STREET ROAD • PO BOX 220 • SKANEATELES FALLS, NY 13153-0220 • PHONE [800] 535-6663 OR [315] 685-4560

Welch Allyn[®]