



Welch Allyn MFI Sōlarc[®] Headlight System

Easy to wear. Hard to wear out.



WelchAllyn[®]

Advancing Frontline Care[™]


Added adjustments for better comfort

Fine tune the two ratchet adjustment systems until the soft, flexible headpiece fits just right.



Spot adjustment
from 20mm to 75mm at
16" operating distance

Luminaire improves your view
Precision optics are securely aligned in a
lightweight alloy luminaire body.

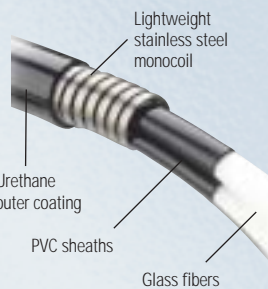

Lightweight stainless steel monocoil
Urethane outer coating
PVC sheaths
Glass fibers

Our triple-shielded fiber optic bundle comes with an ironclad guarantee

We're positive our patented interface and triple-shielding will protect the glass fibers from breaking. So positive, in fact, that we're making this exclusive Ironclad Guarantee. Our MFI fiber optic bundle will maintain 90% of original transmission capacity throughout the first 18 months of use or we'll replace it. Free of charge.



Sleek design
complements any setting

Manual dimmer
provides control of light output
from 0% to 100%

Multiconnection turret for Sölarc light source

The Sölarc light source is also available with a multiconnection turret that lets you plug in most Storz, Olympus, Wolf, and ACMI fiber optic bundles.



Xenon color at a halogen price

Halogen headlights offer yellow, less intense illumination. Sölarc delivers the pure white light needed to make precise tissue color identifications.



Halogen Lamp
(3,200°K)

Welch Allyn
Sölarc Lamp
(5,500°K)

Compact Mobile Stand

Mobile stand securely holds light source in place. It offers storage for your headlight, fiber and spare lamp. The tilt-free design has a compact, 19" diameter footprint.



Welch Allyn MFI Sölarc Headlight

Order Number	Description
49520	MFI Sölarc Headlight System (light source, 8' fiber bundle, headband/luminaire)
49520TW	MFI Sölarc Headlight System with Turret (accepts: Wolf, ACMI, Storz, Olympus)
49525	MFI Sölarc Headlight System, International Power Cord (specify country)
49500	MFI Sölarc Light Source only
49500TW	MFI Sölarc Light Source with Turret (accepts: Wolf, ACMI, Storz, Olympus)

Order Number	Description
49505	Sölarc Light Source only, International Power Cord (specify country)
49540	MFI Headband/Luminaire
49541	Bifurcated Fiber Optic Bundle, 10'
49543	Bifurcated Fiber Optic Bundle, 8'
49550	MFI Sölarc Mobile Stand
09500-U	MFI Sölarc Replacement Lamp



Welch Allyn
Skaneateles Falls, New York 13153-0220
(p) 315.685.4100 • 800.535.6663 (f) 315.685.2174
www.welchallyn.com



Advancing Frontline Care™